

*DCI* edge



DCI EDGE SERIES 5 OPERATORY EQUIPMENT



*CELEBRATING*  
*40 YEARS*  
OF PASSION  
AND PROMISE

DCI is a multigenerational family-run company based in Newberg, Oregon. We believe in managing every aspect of the business around one simple goal: **making our customer's life easier.**

This focus has made DCI Edge the fastest growing dental equipment brand in North America. Our passion is designing and building the best dental equipment in the industry. Our promise is to manufacture with an unwavering focus on reliability and optimal value.





CRAFTED WITH  
DECADES OF  
EXPERIENCE

*DCI* edge



CHAIR CITRUS ULTRALEATHER™ UPHOLSTERY

844-854-8555 | [dciedge.com](http://dciedge.com)








# FINANCIAL FREEDOM

It's your practice. It's your budget. Keep your options open. DCI Edge provides you with all the equipment you need at a price point that allows you to optimize the design and function of your operatory without compromising quality.



## REINVEST

Upgrade your practice with the latest and greatest technology, and secure your edge in the ever-changing tech world.

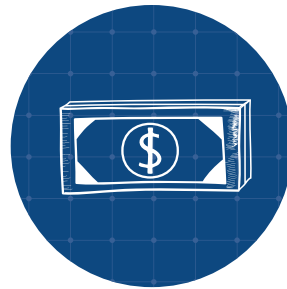

## SAVE

There is peace of mind that comes with extra money in the bank. Save for a rainy day or for your retirement plan.



## CHOICE

Choose where to allocate your funds to make your life easier and keep your practice running smoothly.




# GET COMFORTABLE STAY HEALTHY



## ERGONOMICS

The thin, narrow back design and dual articulating headrest allow for ideal positioning for the provider and patient, ensuring enhanced comfort and efficiency for the entire dental team.

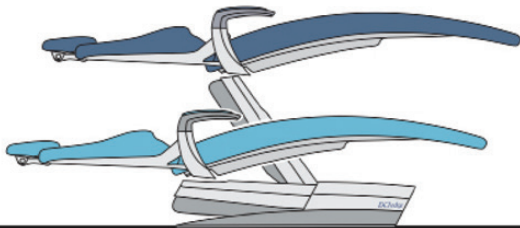

## SIT OR STAND?

The Series 5 chair height adjusts between 16" and 31" so you can comfortably practice both sitting or standing.

## ACCESS TO THE ORAL CAVITY


Unique contoured back design ensures full access between the 3 and 9 o'clock working positions while comfortably supporting the patient.



SERIES 5 CHAIR & STOOL

TIN MAN PEARLIZED ULTRALEATHER  
UPHOLSTERY





## ARE YOU IN THE RIGHT POSITION?

Optimal positioning near the patient ensures a healthy posture, reducing pain and fatigue.



# COMFORT THAT KEEPS PATIENTS SMILING



CHAIR & STOOL APRICOT ULTRALEATHER UPHOLSTERY





**ERGOBACK PLUSH**

## PATIENT SUPPORT

The innovative back cushion and headrest scoop designs provide patients with support regardless of their height.

## SOFT START AND STOP

Adjust your chair without disrupting your patient with smooth start and stop technology.

## DESIGN APPEAL

Available in ergoback or narrow back, and plush or asepsis styles, each with cushioned armrests and a wide variety of color options.

## MULTI-LAYER SYSTEM

Four layers of specialized foam with varying levels of density deliver premium comfort and support for all body types.




**NARROW BACK ASEPSIS**



**NARROW BACK PLUSH**



**ERGOBACK ASEPSIS**

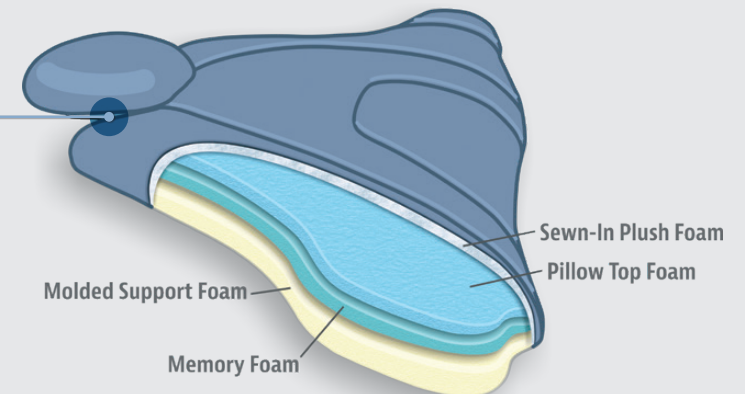

CHAIR

BAYOU ULTRALEATHER UPHOLSTERY



**BEAUTIFUL  
INSIDE AND OUT**








## RELIABLE AND BUILT TO LAST

Our design philosophy blends innovation with high quality construction. We are proud to show you what's under the hood.

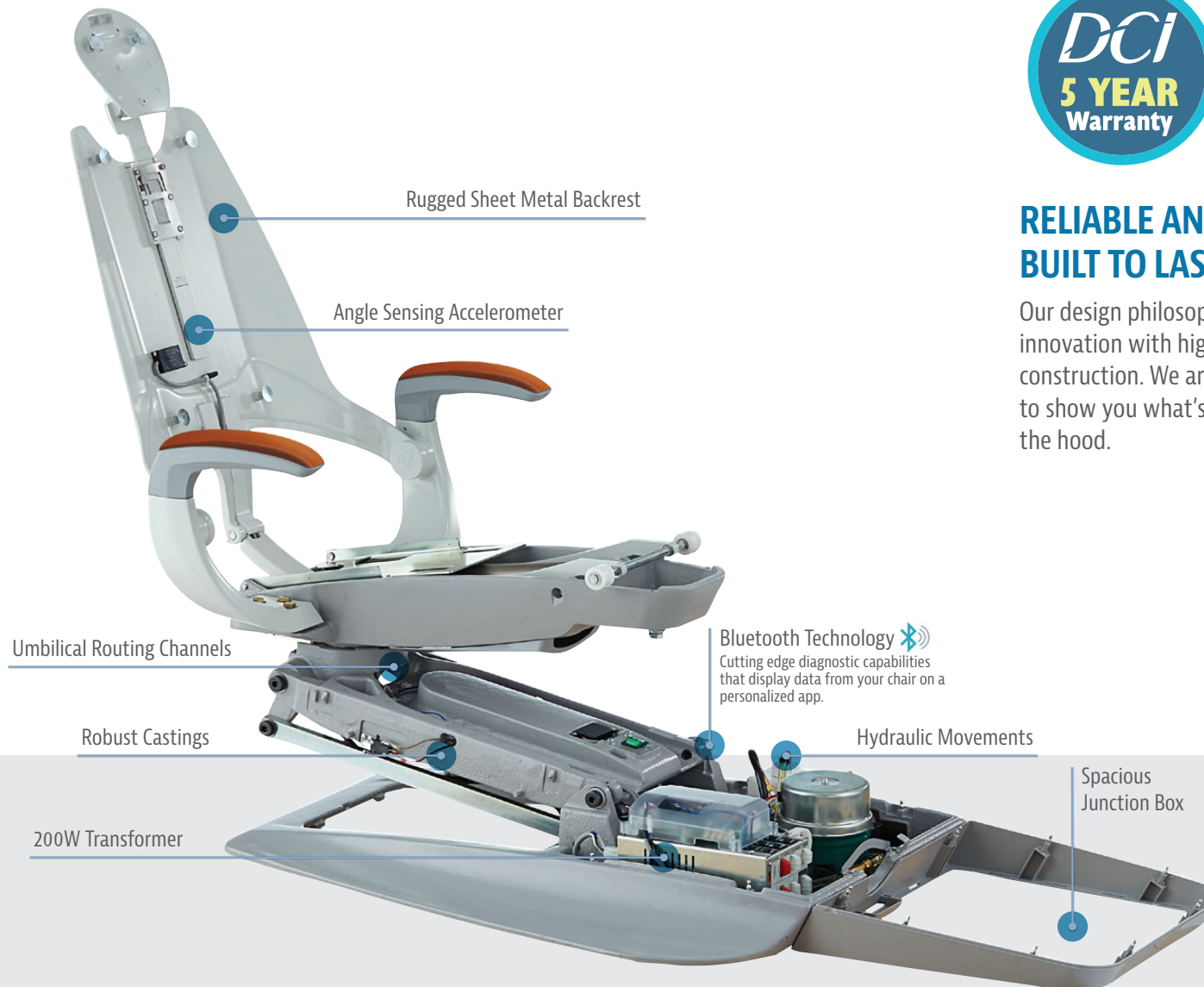



# PERFORMANCE MEETS DESIGN

**NSK**



## EXCLUSIVE EXTENDED WARRANTY

Ask about the DCI Edge exclusive NSK pricing and extended warranty on a variety of air and electric attachment bundles.



## INTEGRATED AND READY TO GO

Easily access essential clinical and diagnostic devices such as the Dentsply Sirona Cavitron system, the ACTEON® SOPRO family of intraoral cameras, piezoelectric scalers, and curing lights.

## DUAL INTEGRATED TOUCHPADS AND BRAKE HANDLES



## ERGONOMIC ACCESS TO INSTRUMENTS

Additional head pivot swivels to improve access and reduce unnecessary movement.



Lights On/Off



Return To Last Position

01


Programmable Buttons for Preset Positions



Backrest Incline and Recline Chair Up and Down

## INTERNAL COMPONENTS


Reliable pneumatic components and ample room for high tech integration. As the largest manufacturer of dental parts in the world, DCI is known for quality parts and components.

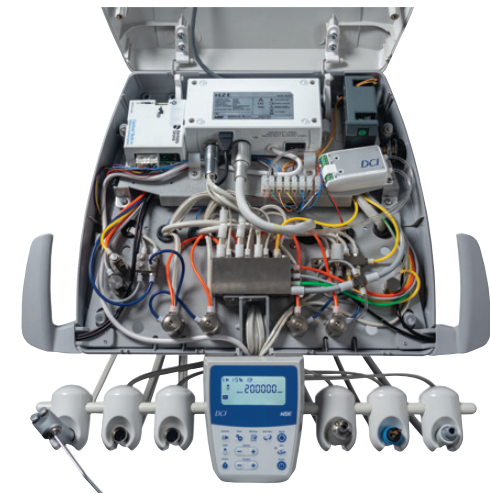




## INTEGRATED TECHNOLOGY

Experience award-winning design featuring WAVEONE™ technology. The NSK NLZE Electric Motor is engineered with superior power, greatly expanding dental treatment possibilities.

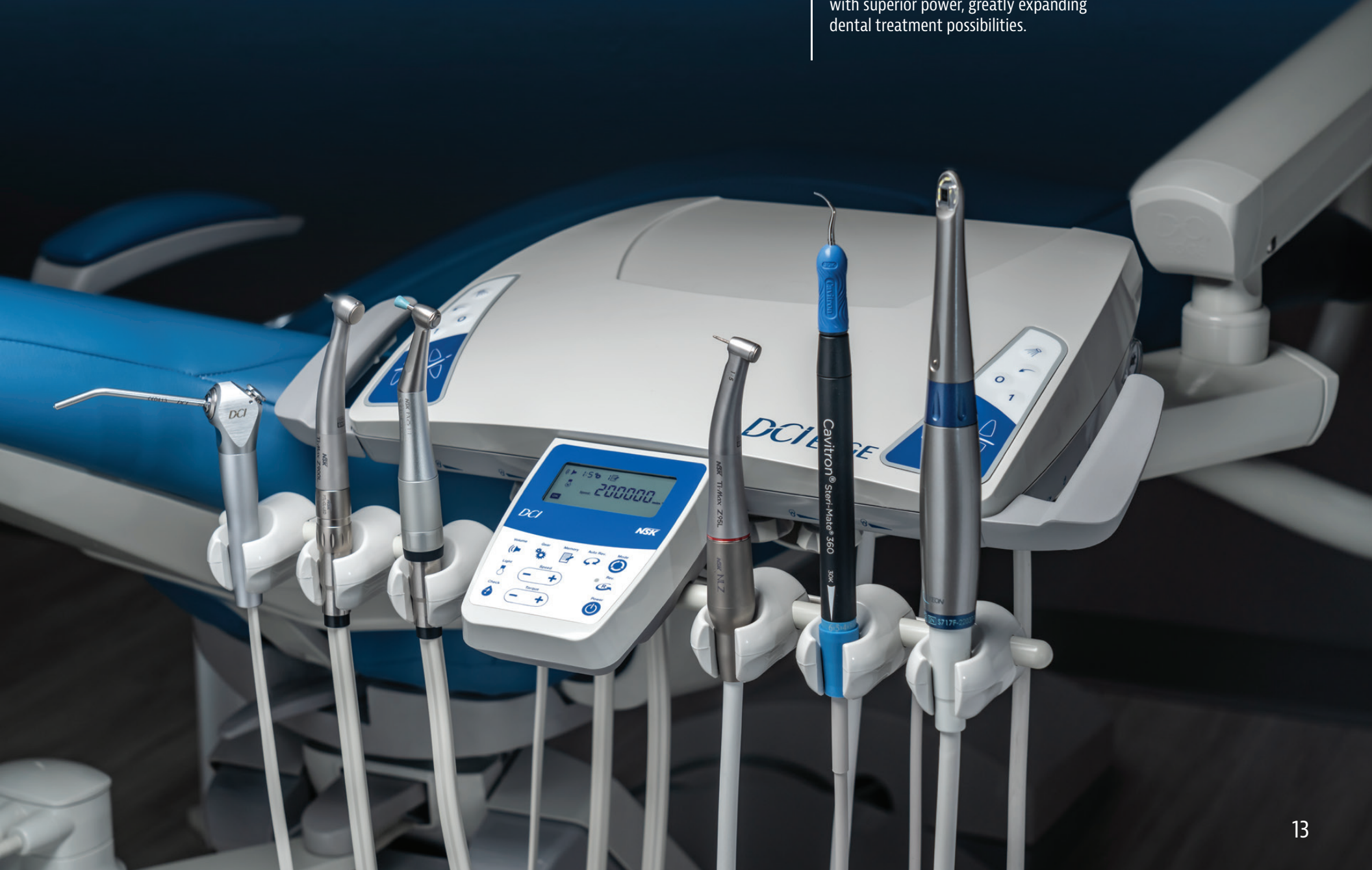



# INNOVATIVE REFLECTIVE TECHNOLOGY INCREASED VISIBILITY

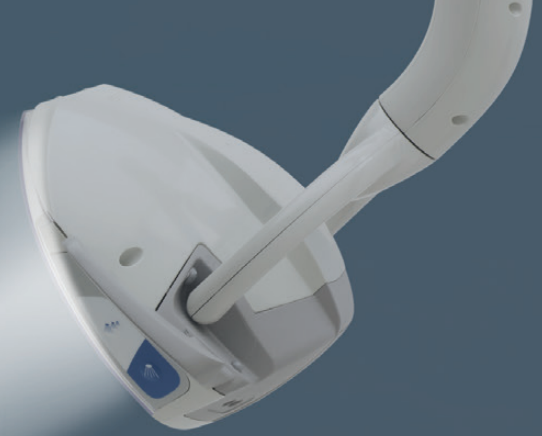

## INCREASED VISIBILITY

DCI delivers cutting edge illumination technology by utilizing over 500 reflective facets to bounce light back to the oral cavity. The result is a uniform light pattern.

## LESS EYE FATIGUE

Optimal light performance reduces shadows, improving illumination of the oral cavity.

## ACCURATE SHADE MATCHING

5000 Kelvin color temperature with a color rendering index (CRI) of 90 produces a neutral white light, which is ideal for color matching at any intensity level.

## NO-CURE SETTING

Filters out UV light to prolong the working time of photo-initiated resins.





## PROGRAMMABLE NO-TOUCH SENSORS

Control the light without contact.



## INTEGRATED TOUCHPADS



No-Cure




Intensity



On/Off

## ACCURATE SHADE MATCHING 5000 K

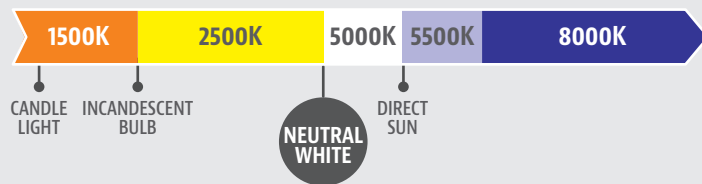

## MOUNTING CONFIGURATIONS



TRACK MOUNT LIGHT



POST MOUNT LIGHT



CEILING MOUNT LIGHT



WALL MOUNT LIGHT  
(Cabinet mount also available)

# SYSTEM SOLUTIONS

## MAXIMUM FUNCTIONALITY



### SWING MOUNT DELIVERY SYSTEM

Delivery system and light easily convert from left- to right-handed operation. The compact swing radius optimizes space, making it perfect for operatories of any size.



SHOWN WITH REAR CUSPIDOR





SHOWN WITH CUSPIDOR

## POST MOUNT DELIVERY SYSTEM

The post mount delivery system features an integrated vacuum trap and water system, along with air and water quick connects.



# SLEEK CABINETRY DESIGN TO COMPLETE YOUR OPERATORY

Designed for optimal organization and efficiency, DCI Edge Cabinetry combines the work and storage space you need with a beautiful modern aesthetic. Purposeful and streamlined construction at a price point that offers the best value for your budget.



CHAIR & STOOL  
SOLID SURFACE  
CABINETS

SKY ULTRALEATHER UPHOLSTERY

MORNING ICE

UPPER : ITALIAN SILVER ASH, LOWER: STERLING ASH









# DCI EDGE CABINETS

## QUALITY, RELIABILITY, AFFORDABILITY

### HIGH QUALITY CONSTRUCTION

Structural steel support with a 6-inch sub-base provides durability while lifting wooden components off the floor and away from moisture. Thermal-fused edge-banding makes doors and drawers water resistant and protects edges from damage.

### PRECISION MANUFACTURING

Designed with cutting-edge software, and cut with state-of-the-art CNC Machinery, DCI Edge cabinetry represents quality at the highest level.

### RELIABLE COMPONENTS

Engineered with top-of-the-line components that stand the test of time.

### FUNCTIONAL AND AFFORDABLE

Built to support all operatory needs, including units and lights, at a competitive price point.

### ATTENTION TO DETAIL

Enjoy detailed features that make routine tasks more efficient: soft-close doors and drawers, undermount sinks that are easily cleaned after each patient, powder-coated steel sub-bases that stand up to repeated cleaning.

### COMPLETE SOLUTIONS


From side cabinets to 12 o'clock workstations, center island consoles to sterilization centers and lab casework, DCI Edge cabinetry offers complete solutions to meet your individual needs.





CHAIR & STOOL  
SOLID SURFACE  
CABINETS

CITRUS ULTRALEATHER UPHOLSTERY  
SILVER SMOKE  
UPPER & LOWER: MAGNOLIA, ACCENT: NIGHTFALL

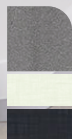



# BEAUTY MEETS FUNCTION

With DCI Edge cabinetry, aesthetics are not compromised for utility. With a designer's attention to visual appeal, our cabinetry provides ergonomic access to materials and instrumentation, ample workspace to accommodate any clinical procedure, and optimal storage space. Additionally, the units are easily adjustable for left- or right-side usage.







**36" 12 O'CLOCK CABINET**  
WITH UNDER COUNTER MOUNT  
DELIVERY SYSTEM



**42" 12 O'CLOCK CABINET**  
WITH CABINET MOUNT  
DELIVERY SYSTEM



**42" 12 O'CLOCK CABINET**  
EXTENDED UPPER WITH REAR  
FLOOR-MOUNTED DELIVERY SYSTEM



# DESIGNED FOR DENTAL ORGANIZATION AND EFFICIENCY

## FLEXIBLE DESIGN

DCI Edge provides modular customization to accommodate even the most challenging spaces and optimize operator organization.

## TECHNOLOGY FRIENDLY

Tech-friendly modules allow for easy and efficient computer and technology integration. Designed with paperless offices in mind.

## SOLUTIONS FOR PRODUCTIVITY

Adjustable pull-out ergo shelves, tip-out bins, flexible monitor mount, integrated electrical and USB ports, and colored drawer organizers ensure that everything you need is efficiently placed within reach.







# INTEGRATED DESIGN YOUR CHOICE OF OPTIONS

Choose the configuration and any additional customization options that allow you and your team to access everything you need without leaving the chair.

## REAR FLOOR-MOUNTED DELIVERY UNIT FEATURES

- Accommodates ancillary integration, now or in the future
- Mobile surface has an adjustable height range of 3"
- Unique design accommodates easy positioning for left- or right-handed practitioners
- 6" side to side movement without encroaching on operator space
- Alignment of collection trap and water bottle enable convenient service position for Assistant
- Assistant touchpad pivots 180 degrees
- Work surface tray holder provides additional space for supplies





## 12 O'CLOCK FREE STANDING CABINET FEATURES

- Assistant's work surface can increase effortlessly to accommodate any clinical procedure
- Colored drawer organizers optimize storage space
- Ergonomic access to tubs, trays, cassettes, materials, and instrumentation
- Upper doors feature textured glass with LED lighting
- Pass-through upper with textured glass inserts
- Mid-section power station with tip out bins. Two access doors with hidden duplex outlets for clean cord management
- Swivel steel powder-coated ergo shelves to hold tubs and trays, or act as a work surface. Robust, die cast arms that hold up to 10 pounds
- Standard 24" or extended 30" upper cabinets
- LED task lighting and clock timer
- Dual glove dispensers
- Removable back panel with freestanding 12 o'clock
- Lower center module with CPU storage, bi-fold door and adjustable shelf
- Lower Assistant's module side opening bi-fold door and pull-out drawers
- Easy-glide monitor mount can be installed at top or bottom of upper cabinet for convenient viewing
- Optional nitrous flow meter cut-outs available
- Accent colors on mid-section



# CABINETRY STREAMLINED FOR STERILIZATION





## FLIP DOWN WRAP STATION

Two conveniently spaced drawers designed to pull out and flip down, ensuring ample space to wrap cassettes of any size.



## SYSTEMATIC DECONTAMINATION

Optimize workflow with specialized waste containers, stainless steel undermount sink, dual integrated drain boards and conveniently placed glove dispensers.

## FACTORY INTEGRATED ACCESSORIES

Our designers and engineers have curated a list of the best options in dental sterilization. Factory installation provides optimal placement and smooth workflow.



## STERILIZATION TOWERS BUILT TO LAST

Built with waterproof materials that guarantee protection from water damage. Multiple tower options accommodate your specific sterilizer selections.

## SIMPLE STERILE ORGANIZATION

Pull-out storage with wire racking holds numerous cassettes, tubs, trays and wrapped instruments. Adjustable shelves allow for fast and easy management.



## SPECIALIZED LIGHTING FOR VISUAL WORKFLOW

Red lights indicate contaminated storage while blue lights illuminate sterile storage. Additional task lighting brightens the entire work surface for maximum visibility.


# CABINETRY STREAMLINED FOR WORKFLOW



RECEIVING, CLEANING  
& DECONTAMINATION



PREPARATION  
& PACKAGING

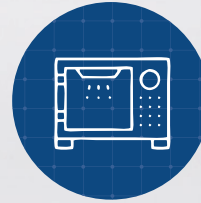

STERILIZATION



STERILE STORAGE  
& MAKE READY





## A WORKSPACE BUILT JUST FOR YOU

Sterilization is one of the most important steps in a practice, and one size never fits all. With DCI Edge, you can create a Sterilization Center that fits your unique needs. All options feature a process-based flow that keeps productivity high and risk low.

## STERILIZATION WITH STANDARDS

Designed to follow CDC guidelines with exceptional focus on infection control and patient safety.


## STERILIZE WITH STYLE

Fully functional clock timer enhances productivity, while stainless steel valance adds to the beautiful design aesthetic.


### STRAIGHT RUN




### GALLEY STYLE

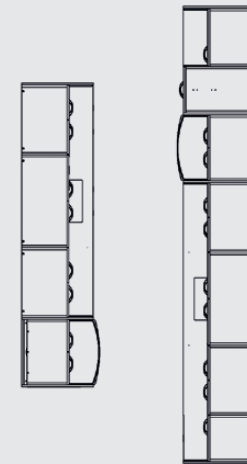

### L SHAPE

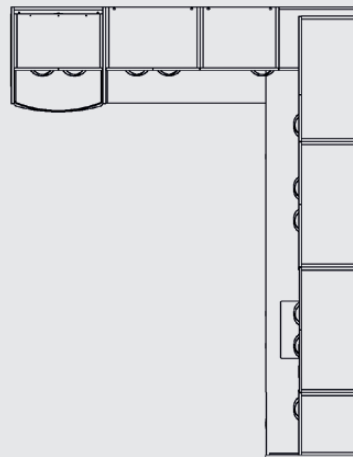

### U SHAPE

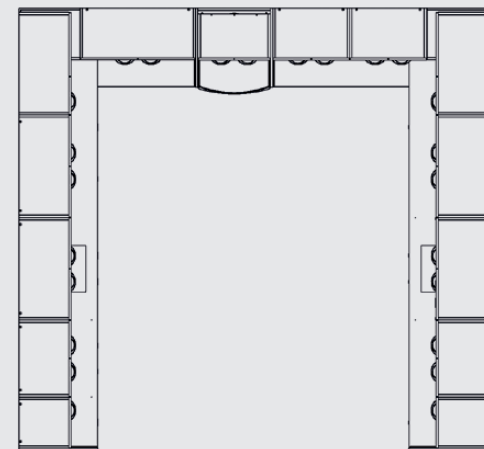



## EXPERIENCE DCI EDGE

3 SHOWROOMS

2 MANUFACTURING FACTORIES



**DCI Edge Showroom &  
Manufacturing Facility  
Oregon, USA**

120 S Elliott Road  
Newberg, OR 97132



**DCI Edge Showroom  
and NSK America  
Corporate  
Headquarters**

1800 Global Parkway  
Hoffman Estates, IL 60192



**DCI Edge Showroom &  
Manufacturing Facility  
North Carolina, USA**

9201 Stockport Place  
Charlotte, NC 28273

**dciedge.com**

**equipment@dcionline.com**

**1-844-854-8555**



50169/August2023-003  
©2023 DCI

