
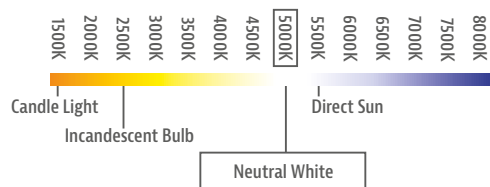


DCI *edge*

Dental Lights

Cutting-Edge Illumination Technology

DCI Edge Series 4 and Series 5 lights feature breakthrough LED technology with over 500 reflective facets to eliminate shadowing and deliver next-level lighting performance.

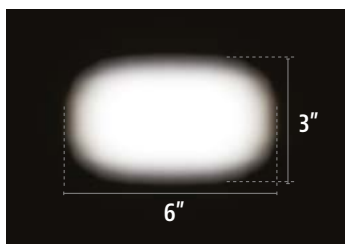

Reflective LED Technology

DCI delivers cutting-edge illumination through 500+ reflective facets that direct light precisely into the oral cavity, creating a uniform 3" x 6" light pattern that minimizes shadows and maximizes visibility.

Accurate Shade Matching

A 5000 Kelvin color temperature with a color rendering index (CRI) of 90 produces a neutral white light that ensures accurate color matching at any intensity level.

3rd Axis Rotation

Improves light access in all quadrants of the oral cavity.

Durable Design

Die-cast construction provides long-lasting durability.

No-Touch Sensors

Control the light without contact using programmable no-touch sensors, available in the Series 5 LED light.

Compare DCI Edge Lights

FEATURES	SERIES 5	SERIES 4
LED Technology	✓	✓
Reflective Technology	✓	✓
3rd Axis Head Rotation	✓	✓
3 Intensity Settings (35,000 lux, 26,000 lux, and 17,500 lux)	✓	✓
27" Focal Length	✓	✓
Dual Touch Controls	✓	✓
Offered in Gray or White	✓	✓
3" x 6" Light Pattern	✓	✓
5000K Color Temperature	✓	✓
No-Cure Setting	✓	
Dual No-Touch Sensors	✓	
Programmable Sensors to Control On/Off & No-Cure	✓	

LED Lamp: 24-Volt, 10 watts; Input Power: 115V/230V 50/60 Hz; Power Consumption: 0.2 amp at 115 volts



Configurations in Post Mount, Swing Mount, Ceiling Mount, Track Mount, Cabinet Mount, and Wall Mount styles.



Colors available in both gray and white to match your existing operatory equipment.

Experience DCI Edge Lights in 3D

Step inside four fully equipped operatories with the DCI Edge Visualization Tool and see how the lights look in a clinical setting.

